



NICO myriad-LXII® UI GUIDELINES

September 2022

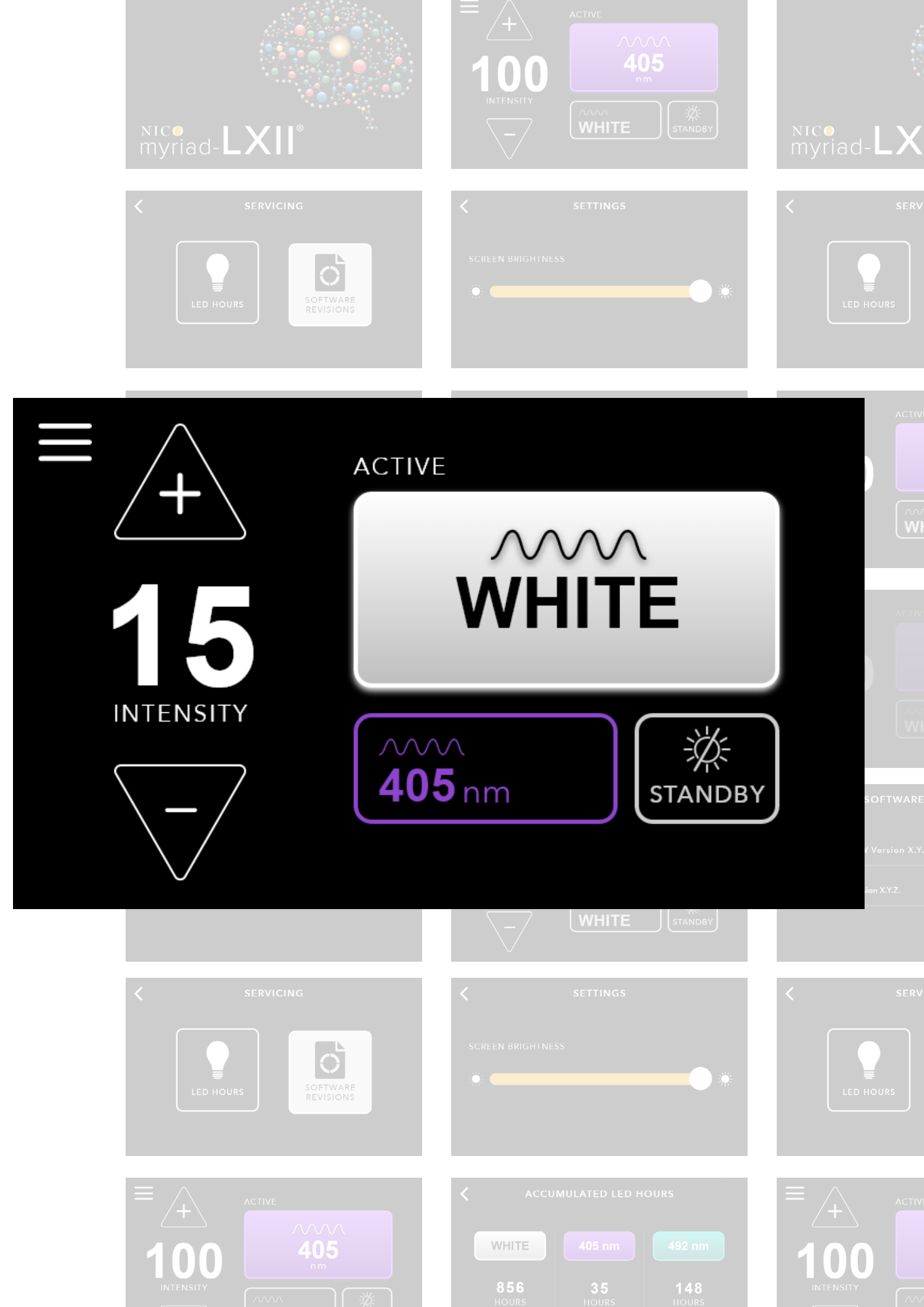

CONTENTS

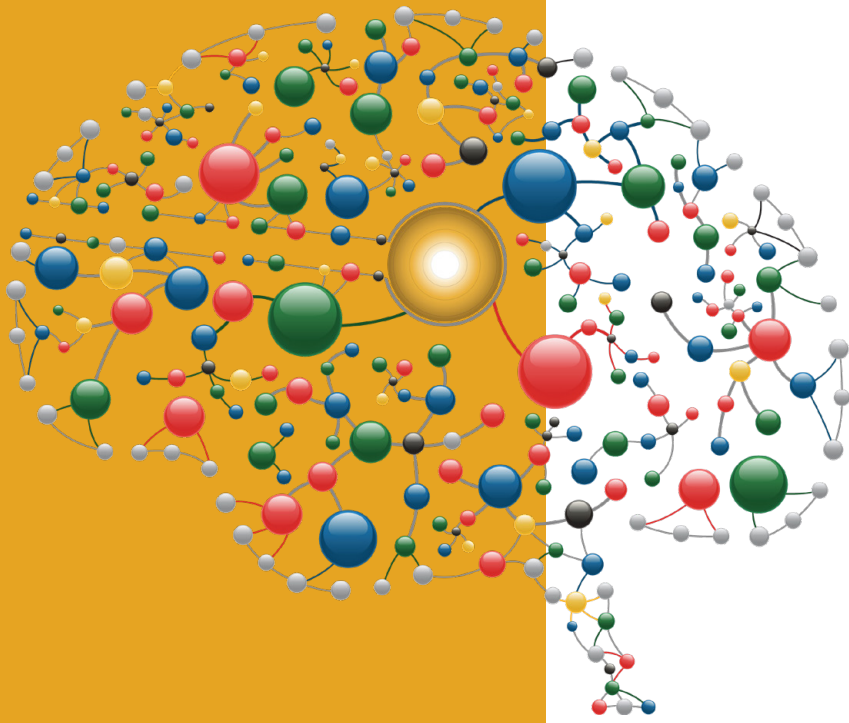

01 INFORMATION ARCHITECTURE

02 USER WORKFLOW

03 TOUCH SCREEN INTERFACE

- Screen Overviews
- Time Delays
- Screen Technology
- Grids + Spacing

04 STYLE IDENTIFIERS

- Font +Typography
- Color + Gradients
- Iconography
- Buttons + Controls

05 PROTOTYPE


06 HARDWARE

01

INFORMATION ARCHITECTURE

01 INFORMATION ARCHITECTURE

The information architecture is a high level diagram which conveys how information is organized within the interface.

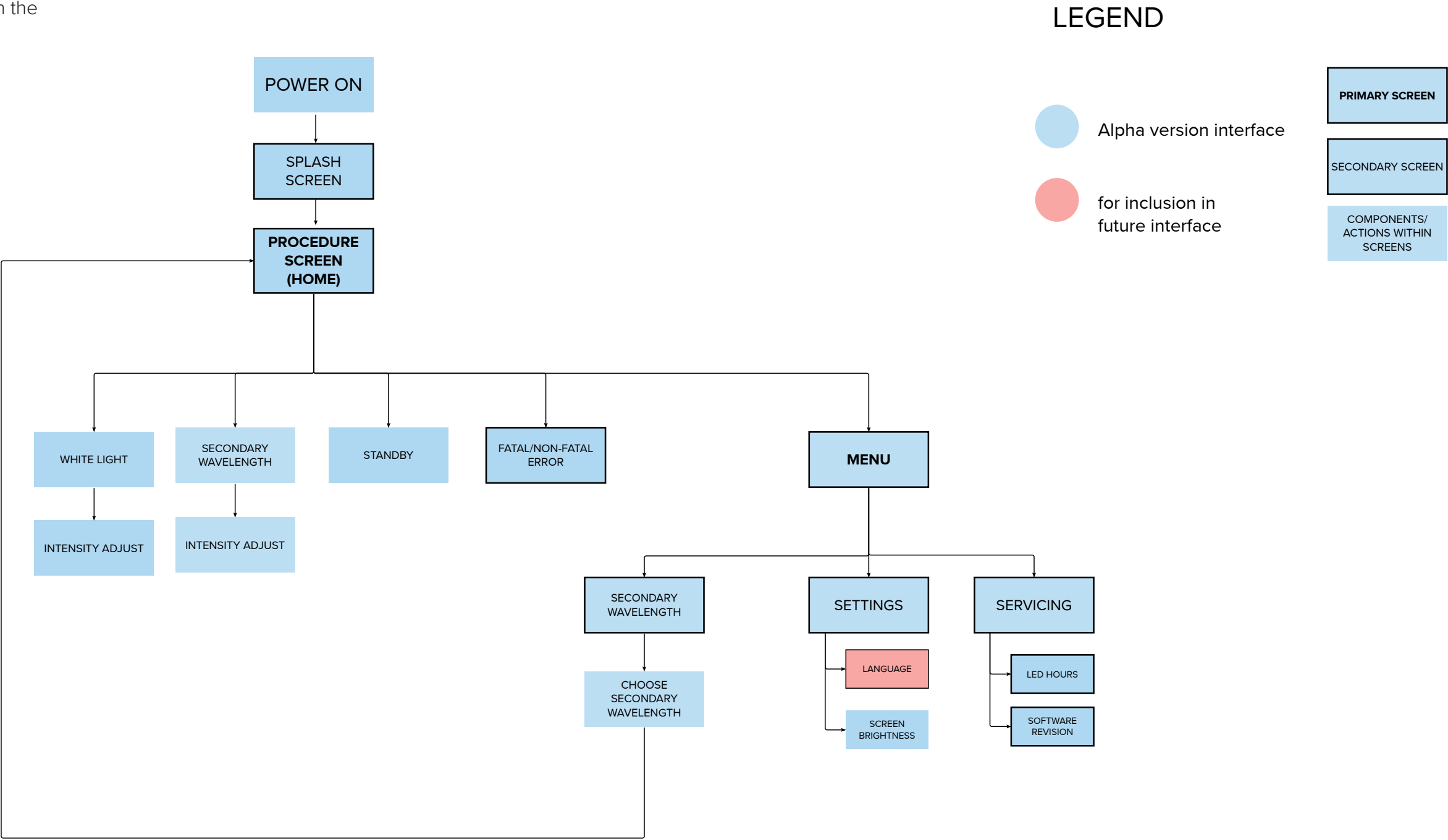


02

USER

WORKFLOW

02 USER WORKFLOW

The user workflow provides an overview of how to navigate through the interface through numerous decision points.

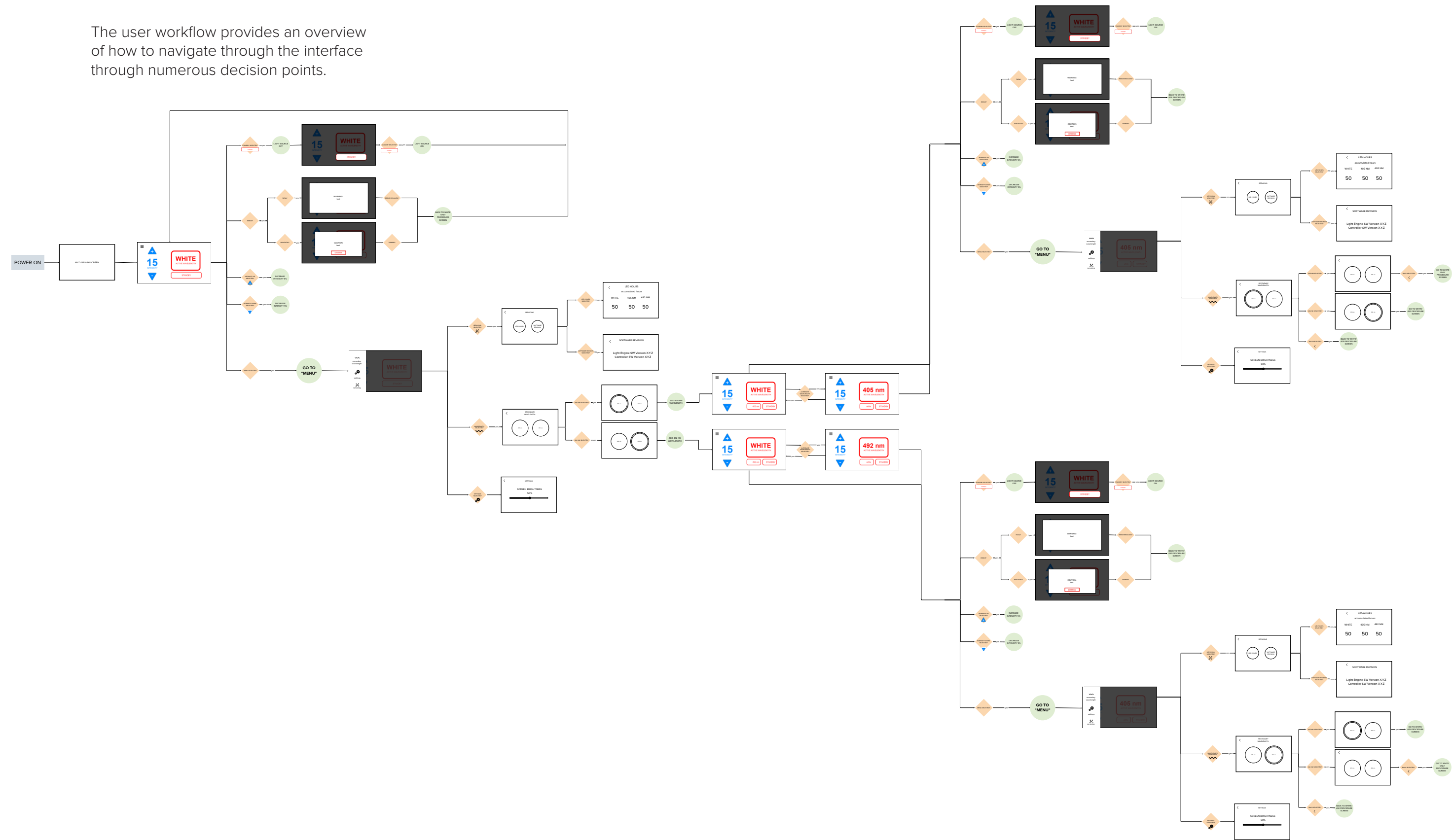

03


TOUCH SCREEN INTERFACE

03 TOUCH SCREEN INTERFACE


PROCEDURE SCREENS

Procedure screens are considered the primary screens in the interface and will be interacted with the most. The procedure screen gives access to intensity adjustment, change of light type, standby mode, and menu options.


WHITE LIGHT ONLY PROCEDURE

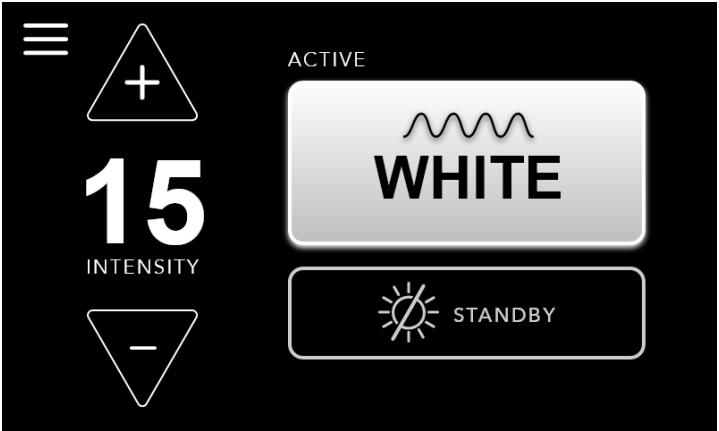

WHITE LIGHT + 405 NM PROCEDURE

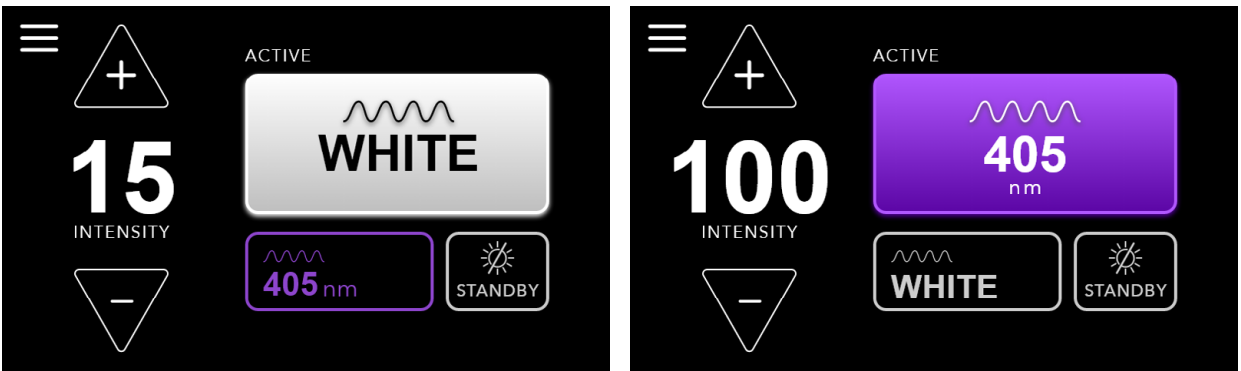

WHITE LIGHT + 492 NM PROCEDURE

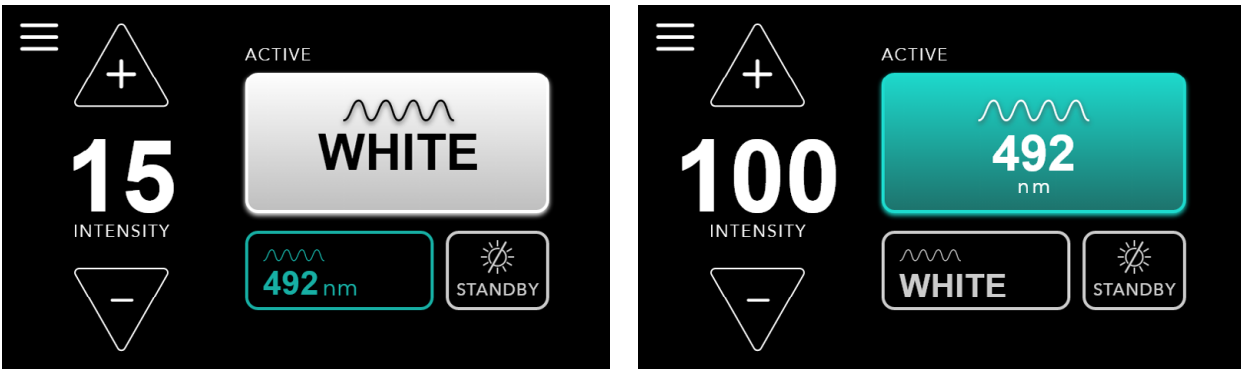

03 TOUCH SCREEN INTERFACE

MENU

The menu is accessed through the hamburger menu icon on the upper left corner of the screen.

The menu is simple, intuitive, and unobtrusive. It slides out from the left side and covers a third of the screen so that wavelength information is still visible. The menu provides access to adding/removing a secondary wavelength, settings (brightness adjustment), and servicing.

Menu icon




SECONDARY
WAVELENGTH


SETTINGS


SERVICING



Menu slides left


03 TOUCH SCREEN INTERFACE

WARNINGS

Caution and warning messages appear automatically when an internal error or problematic temperature increase occurs.

Caution Yellow indicates a cautionary message that is dismissable, while Warning Red indicates a fatal error that must be resolved before the user proceeds.

CAUTION (NON-FATAL ERROR)

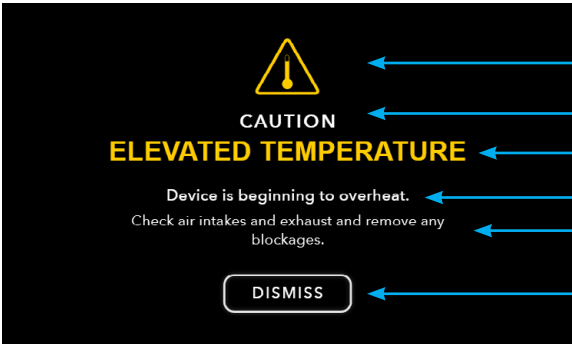

The dialog box is black with yellow text and icons. It features a yellow warning triangle icon at the top. Below it, the word 'CAUTION' is in white, and 'ELEVATED TEMPERATURE' is in yellow. The description text is white, and a 'DISMISS' button is at the bottom.

Icon // indicates type of error

Severity // Caution or Warning

Type // Error Title

Caption // Overview of Error


Description // Details + Instructions

Button // Dismiss error (Caution only)

WARNING (FATAL ERROR)



The dialog box is black with red text and icons. It features a red warning triangle icon at the top. Below it, the word 'WARNING' is in white, and 'TEMPERATURE SHUTDOWN' is in red. The description text is white.

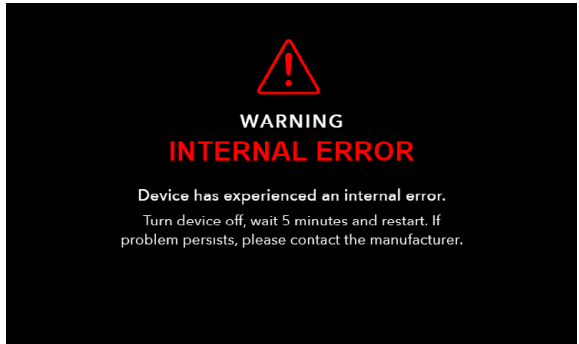


The dialog box is black with red text and icons. It features a red warning triangle icon at the top. Below it, the word 'WARNING' is in white, and 'INTERNAL ERROR' is in red. The description text is white.

03 TOUCH SCREEN INTERFACE


TIME DELAYS

Time delays are incorporated into the following three scenarios to emphasize specific moments and/or interactions within the interface.


STARTUP (SPLASH SCREEN)

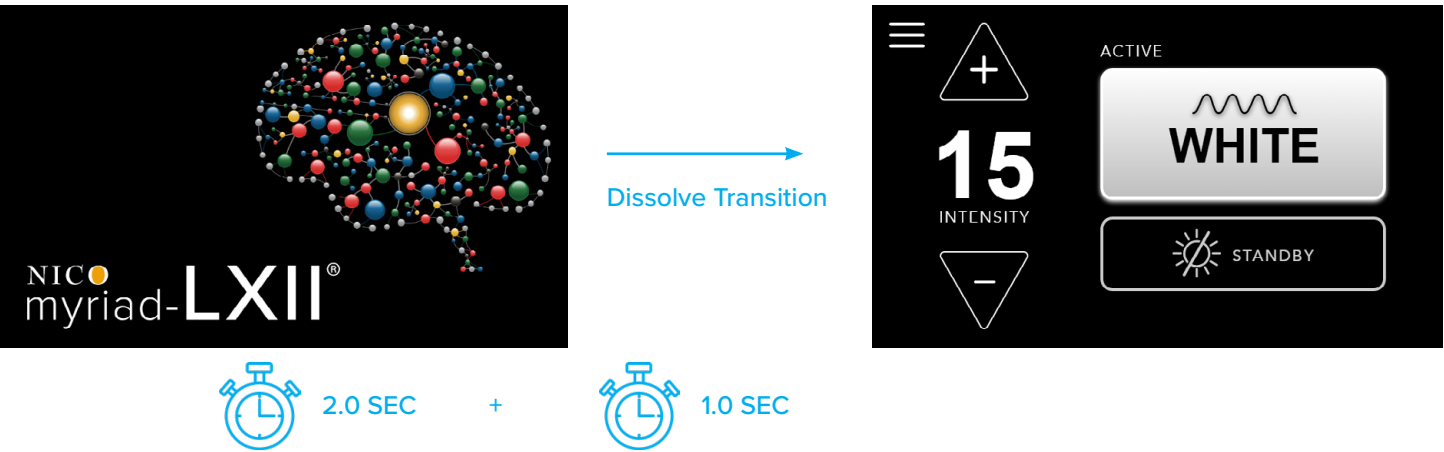

WAVELENGTH SELECTION

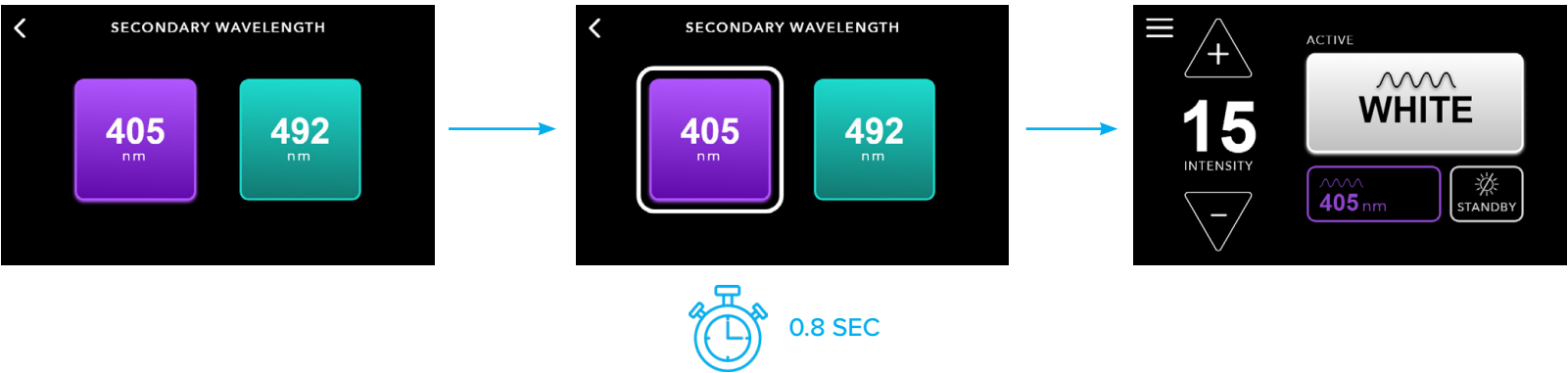

BUTTON STATES

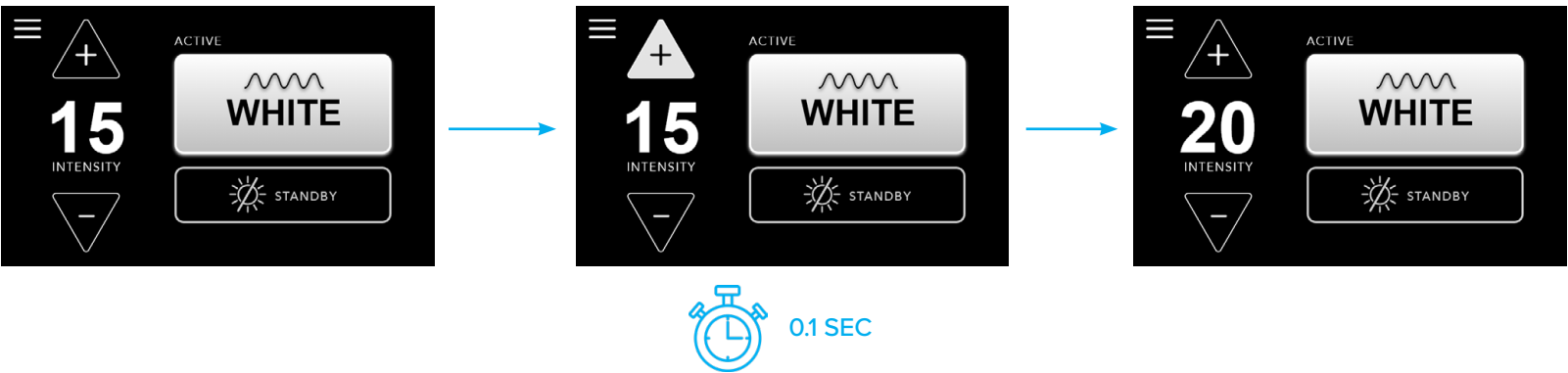

04

STYLE

IDENTIFIERS

04 STYLE IDENTIFIERS

FONTS

The primary font for the Nico Myriad LX II interface is Avenir Next. Avenir Next is a clean, contemporary sans-serif font which is well suited for the operating room environment.

The secondary font is Arial Bold. Arial Bold is a strong, clear and familiar typeface.

ABCDEFGHIJKLMNOPQRSTUVWXYZ
 abcdefghijklmnopqrstuvwxyz
 0123456789 .,:;?&{ }()*%@

Avenir Next, Regular

ABCDEFGHIJKLMNOPQRSTUVWXYZ
 abcdefghijklmnopqrstuvwxyz
 0123456789 .,:;?&{ }()*%@

Avenir Next, Medium

ABCDEFGHIJKLMNOPQRSTUVWXYZ
 abcdefghijklmnopqrstuvwxyz
 0123456789 .,:;?&{ }()*%@

Avenir Next, Demi Bold

ABCDEFGHIJKLMNOPQRSTUVWXYZ
 abcdefghijklmnopqrstuvwxyz
 0123456789 .,:;?&{ }()*%@

Arial, Bold

04 STYLE IDENTIFIERS

TYPOGRAPHY

The Arial font is used primarily for numerical information, including wavelength and intensity value.

The Avenir Next font is used in a variety of weights and capitalizations to create a hierarchy of information on each screen.

01

Arial, Bold, 131 pt

02

Arial, Bold, 66 pt

03

Arial, Bold, 40 pt

TITLE

Avenir Next, Demi Bold, 26 pt, Capitalized

SUBTITLE

Avenir Next, Regular, 22 pt, Capitalized

BUTTON

Avenir Next, Medium 22 pt, Capitalized

Caption

Avenir Next, Medium, 22 pt

Body

Avenir Next, Regular, 20 pt

04 STYLE IDENTIFIERS

COLOR

All screens should always use a solid black background color. Bright white should be used for the majority of text and graphics to ensure the highest contrast possible for optimized visibility. Soft gray should be used in accent instances such as pressed button states. Nico Gold is implemented as an accent color in instances that require more attention and/or interaction.

Wavelength representation colors should be used in any instances in which the 405 nm or 492 nm wavelengths are referenced in order to reinforce the selection of each color.

Caution Yellow and Warning Red should be used during error messages to indicate a non-fatal error or fatal error, respectively.

PRIMARY


screen background



Black
#000000

SECONDARY

text, buttons + accents

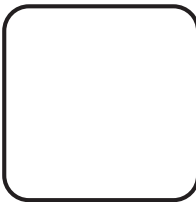

Bright White
#FFFFFF



Soft Gray
#E5E5E5



Nico Gold
#E6A414

TERTIARY

wavelength representations + warnings



Purple
#B157FF
(405 nm)



Teal
#1DDDBCE
(492 nm)



Warning Red
#FF0000



Caution Yellow
#FFCC00

04 STYLE IDENTIFIERS


GRADIENTS

Linear gradations are used to bring visual interest and dimensionality to screens. In this specific interface, gradations fill larger areas of the screen that are used to represent which wavelength is active at a given moment. A filled gradation box behind the name of a specific wavelength indicates that it is active.

Gradations should be vertically linear and move from a lighter tone down to a darker tone.

All gradation boxes should maintain a 4 px outline using the *lightest color* in the gradation scheme.

DIRECTION

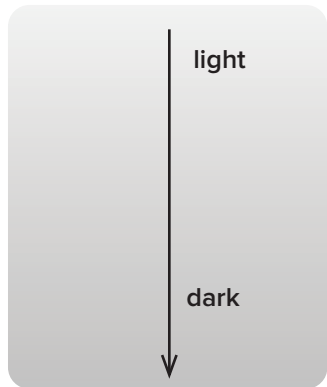

Direction of Linear Gradient

COLOR OPTIONS

White Light



Gradation
#FFFFFF to #BABABA

405 nm



Gradation
#B157FF to #55009F

492 nm



Gradation
#1DDBCE to #1E6B65

OUTLINES




04 STYLE IDENTIFIERS

ICONOGRAPHY

Icons are used to reinforce messages within the user interface. The chosen icons are intended to connect with the user in order to create a more efficient and memorable user experience.

Icons should be as simple as possible to clearly convey the object that they represent. The majority of icons in this set consist of solid filled imagery. In some cases, however, the outlined icon style is used in instances where a solid shape would detract from the overall word or message being conveyed (i.e. warning icons and standby icon).

COMPLETE ICON SET

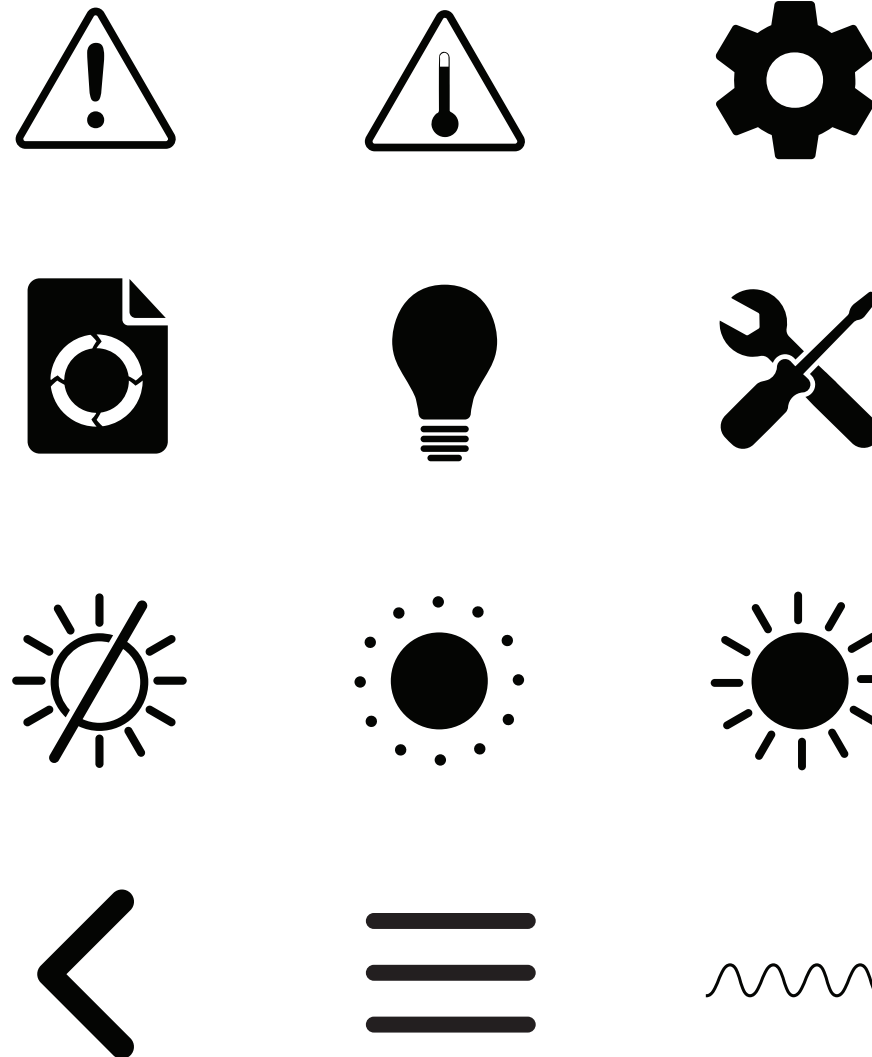


04 STYLE IDENTIFIERS


BUTTON STYLE

Interactive buttons are denoted by a rounded outline shape around a combination of text and/or iconography. Intensity adjustment buttons follow a rounded triangular shape, while all other buttons follow a rounded rectangular shape.

SIZING

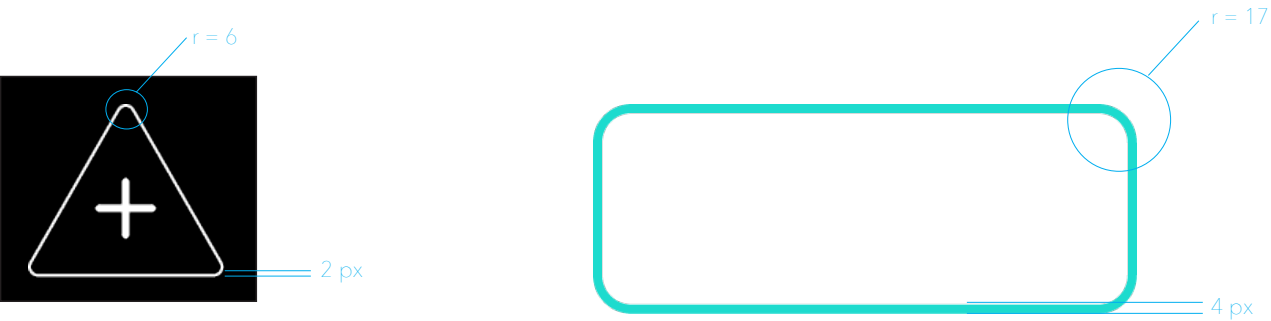

BUTTON STYLE EXAMPLES

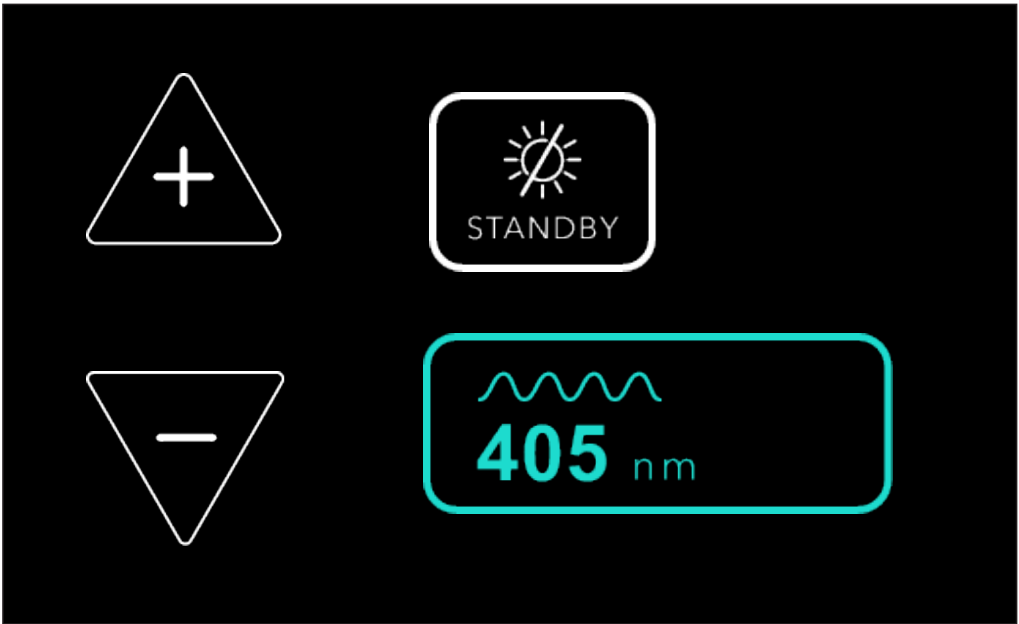


04 STYLE IDENTIFIERS

BUTTON STATES

The standard button style in its default state follows an outlined shape with corresponding inner text and icon color. When a button is pressed, the activated state is a solid fill with Soft Gray, and all inner text and graphics are Black. This instant change provides positive feedback to the user that the specific button has been successfully pressed.

STANDARD BUTTON STATES

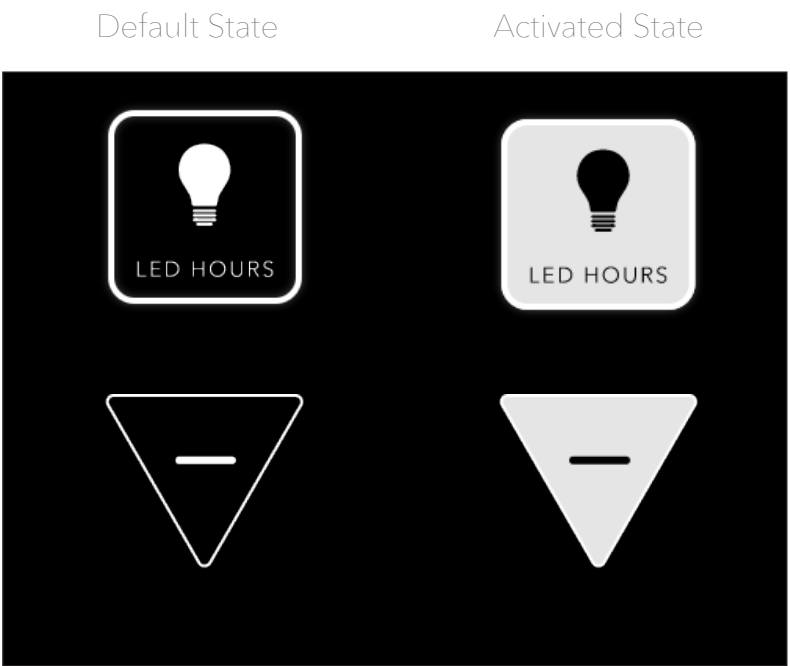

04 STYLE IDENTIFIERS

BUTTON STATES

There are two special cases in which the button states don't follow standard activated procedure.

The standby button, which temporarily turns off the light source when pressed, uses a blinking pattern with Nico Gold in its activated state. This motion and color change is intended to call additional attention to the standby button in order to clearly communicate that standby mode is active. To turn the light source back on, the user presses the blinking standby button again to return to default state.

During wavelength selection, a white outline is used around a wavelength box to communicate that the specific secondary wavelength is selected. This box emphasizes the user's wavelength choice before moving into the procedure screen with the chosen wavelength option.

SPECIAL CASE: *STANDBY MODE*

Default State



Activated State // Blinking Solid and Outline



Default State



Note: All other screen elements, excluding menu icon, are reduced to 30% opacity and remain inactive during standby mode.

SPECIAL CASE: *WAVELENGTH SELECTION*

Default State

Selected Wavelength




04 STYLE IDENTIFIERS

OTHER CONTROLS

Sliders should always display as Nico gold on the activated portion and black, gray or white on the unactivated portion. A white circle should be used as the touch point to interact with the slider.

BRIGHTNESS SLIDER

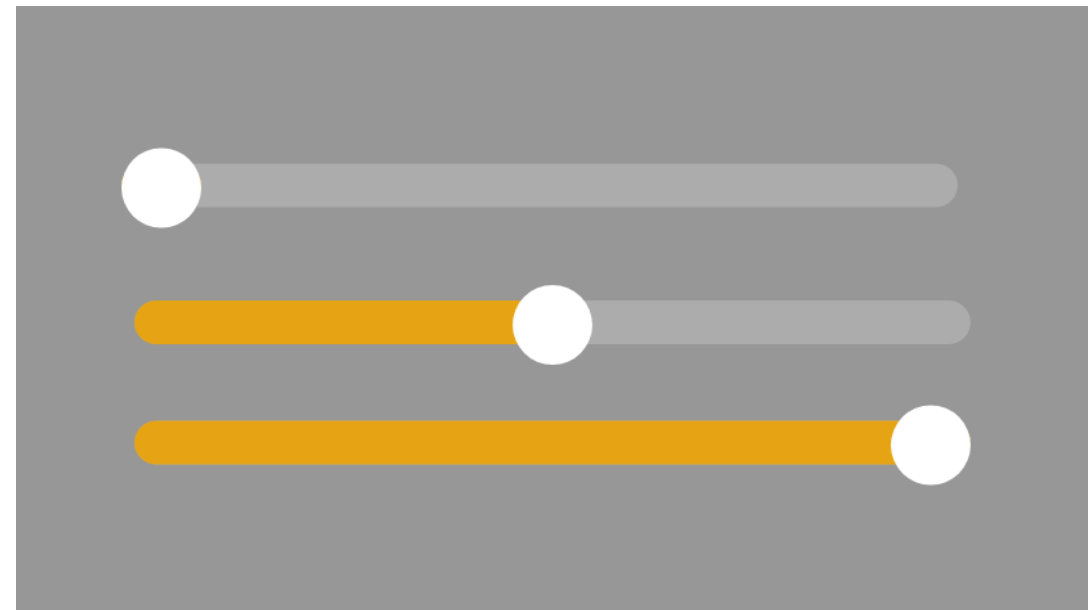

05

PROTOTYPE

05 TOUCH SCREEN INTERFACE

PROTOTYPE

The following link provides access to an interactive prototype for user testing of the interface.

Click Link Below for Prototype:

<https://xd.adobe.com/view/4bc218f4-adee-47ff-87fc-3c5046f6b7af-1e91/?fullscreen&hints=off>


06

HARDWARE

06 HARDWARE

GRIDS AND SPACING

Screen should maintain a minimum margin of 24 px around all four edges to ensure successful touch interaction with the display.

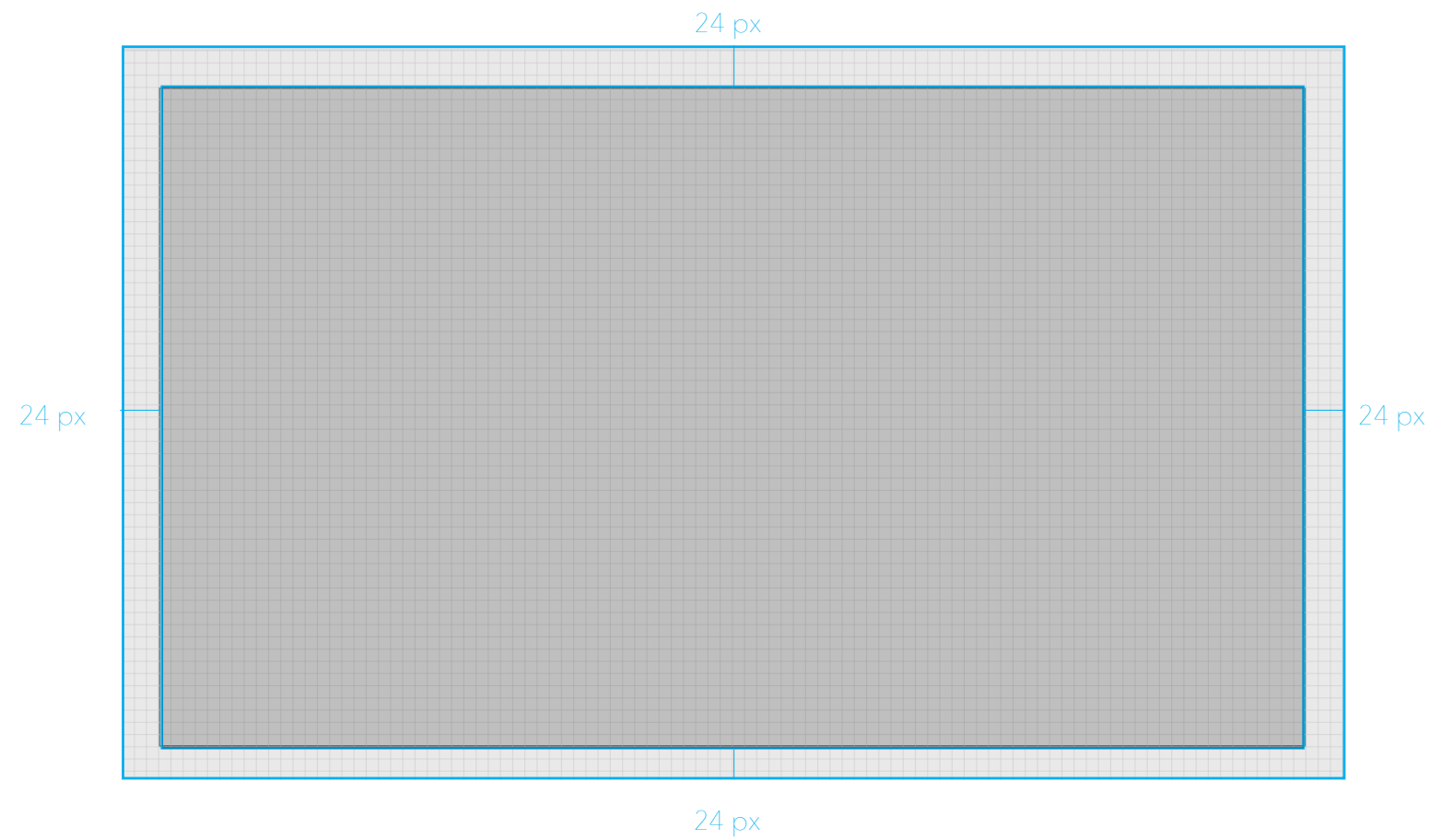


06 HARDWARE

SCREEN TECHNOLOGY

The display for the Nico Myriad LX II is a 7” capacitive touch graphic LCD display screen.

Full Spec Sheet: <https://www.digikey.com/en/products/detail/newhaven-display-intl/NHD-7-0-800480EF-ASXN-CTP/6026270>

SCREEN SPECIFICATIONS

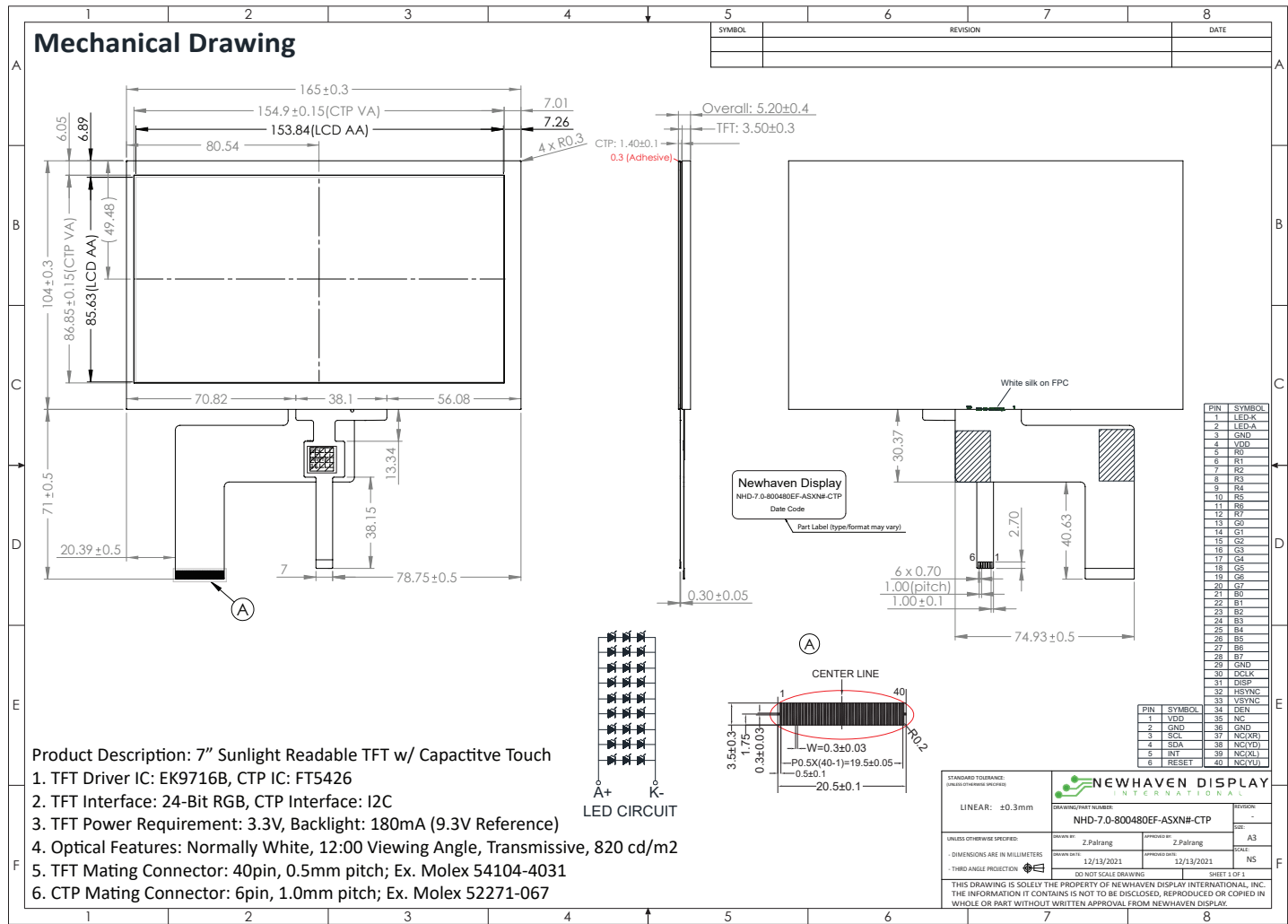

THANK YOU.